


Product Specification

Code	Proximal Diameter(D1) (mm)	Distal Diameter(D2) (mm)	Total Covered Length(L1) (mm)	Delivery System (Fr.)
AUI2412B160	24	12	160	21
AUI2414B140	24	14	140	21
AUI2612B160	26	12	160	21
AUI2614B140	26	14	140	21
AUI2812B160	28	12	160	21
AUI2814B140	28	14	140	21

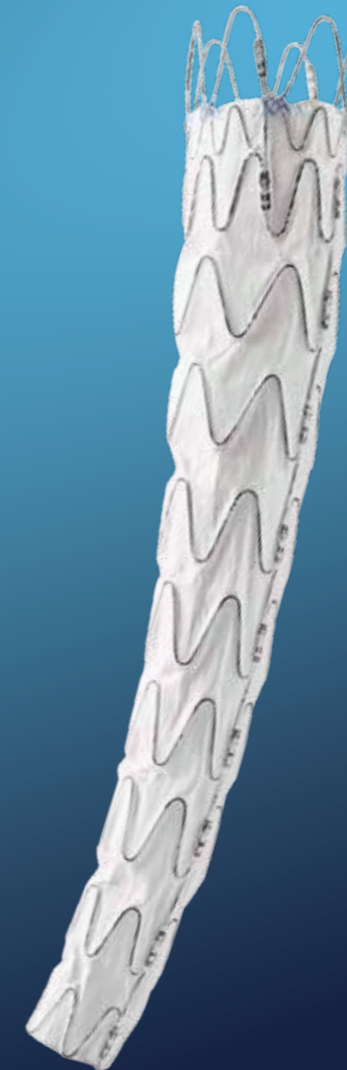


Lifetech Scientific (Shenzhen) Co., Ltd.

Add: Cybio Electronic Building, Langshan 2nd Street, North Area of High-tech Park, Nanshan District, Shenzhen 518057, P. R. China
 Tel: +86-755-86026250 Fax: +86-755-86026251
 Email: marketing@lifetechmed.com
 Web: www.lifetechmed.com

EU authorized representative:
 LifeTech Scientific (Europe) Coöperatief U.A.
 Kruisdonk 64, 6222 PH Maastricht, the Netherlands

2017-EN-A-1.0

Stable support, Precise Apposition

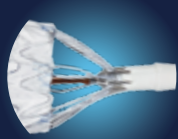

Tip capture mechanism

- provides precise placement and controlled deployment

Proximal mini wave

- enhances apposition to the aortic wall

Connecting bar

- offers axial support
- avoids stent shortening

No suture on the main body

- avoids pinhole leakage (Type IV endoleak)

e-PTFE dual membrane

- excellent biocompatibility, no Type III endoleak



Ankura™ AUI Stent Graft Cases

